


## Connectivity Diagram for Typical Practice

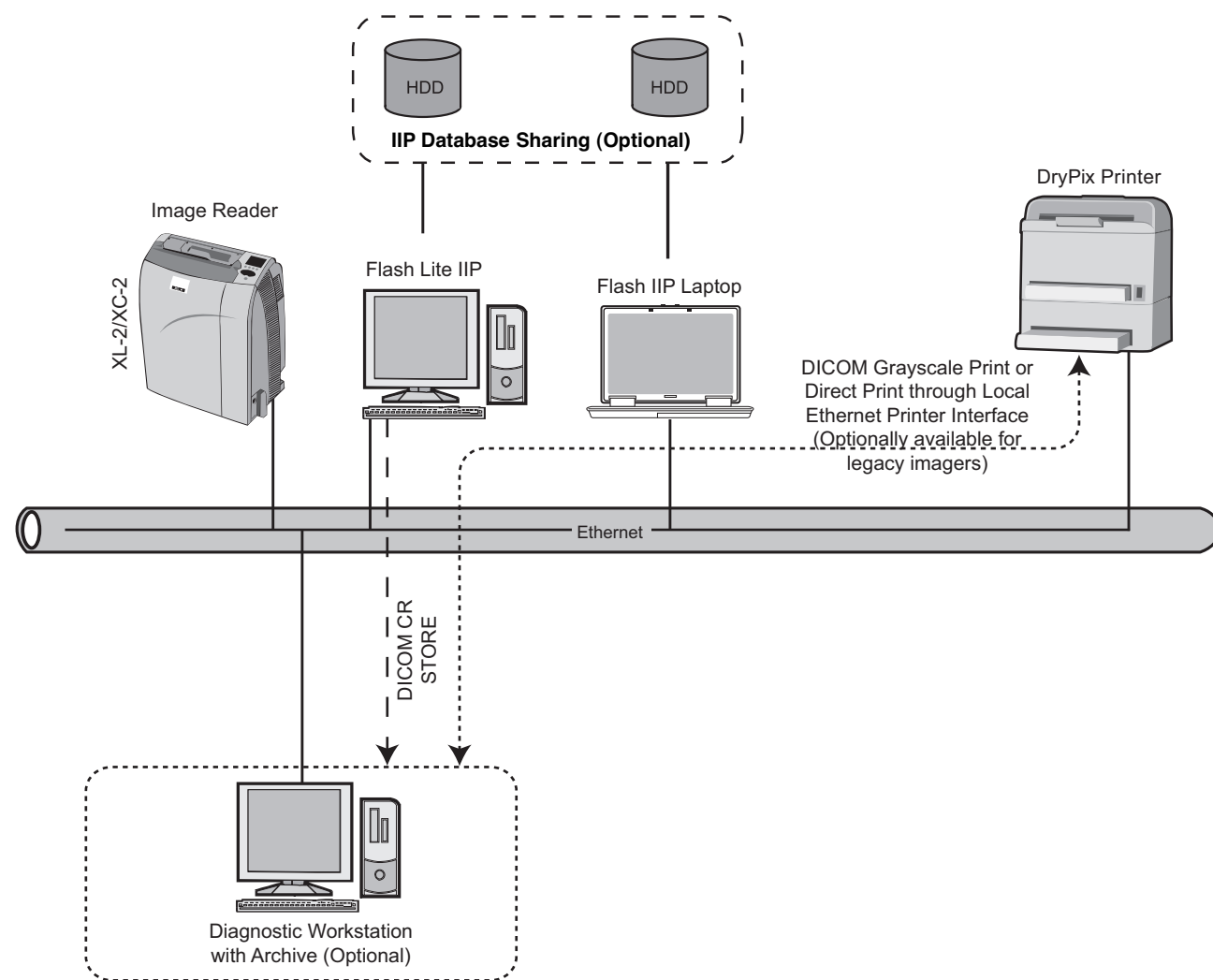

### Notes:

- Network Interface - 100 Base-T, half or full recommended.
- Interface to other FCR Readers and/or Flash IIPs on the network.
- Any FCR Reader can be utilized after registering the patient at a Flash IIP connected to the same network.
- Database Sharing (option) allows sharing exam and image data between Flash IIP Consoles within Database Sharing cluster.
- Interfacing with Legacy FCR equipment - devices that communicate through DMS protocol can be connected to the network printer by adding an optional Fujifilm DryPix Station to convert the DMS to DICOM.

*Specifications subject to change without notice.*

### FUJIFILM Medical Systems USA, Inc.

#### Corporate Headquarters

419 West Avenue  
Stamford, CT 06902-6348  
203-324-2000  
800-431-1850

1055 Stevenson Court  
Roselle, IL 60172-2300  
630-582-2202  
800-323-2546

29012 N. Hancock Parkway  
Building 7  
Valencia, CA 91355-1007  
866-533-FUJI (3854)

www.fujimed.com

### Authorized Fujifilm Dealer:



rev. 2/08

XBUSCH009

# FUJIFILM MEDICAL SYSTEMS



## Product Data

# FCR XL-2 & FCR XC-2

## Digital X-ray with Exceptional Image Quality for Private Practitioners in an Efficient, Compact Package



### Application

The Fujifilm Computed Radiography (FCR) XL-2 and XC-2 are digital radiography systems utilizing storage phosphor Imaging Plates (IP) as an x-ray detection device.

FCR XL-2 and XC-2 units scan and process x-ray image information from an IP that has been exposed using a cassette-type x-ray stand. This system can also print out hard copies on film via a connected image recorder and can transmit digital images to PACS.

The FCR XL-2 and XC-2 provides general x-ray imaging for Private Practitioners. FCR XL-2 and XC-2 systems also include the capability to perform Pantomographic exams.

### Features

Due to the benefits of digital image processing, the system produces radiographs that have a high diagnostic value and are easy to read. Fast Scan (200 $\mu$ ), Standard High Resolution (100 $\mu$ ) and High Resolution (50 $\mu$ ) size using HR IP's is possible with the FCR XL-2.

Reading modes are user selectable. Fast Scan reading provides maximum throughput for the XL-2 for both 14x14" and 14x17" sizes, while standard HQ reading mode provides 10 pixel/mm resolution image acquisition and display for all intake sizes without impacting system speed.

High Resolution (50 $\mu$ ) Scanning is available with the FCR XL-2 for 18x24cm and 24x30cm sized (HR IP's) making it ideal for orthopedic and extremity applications, where seeing fine detail is critical.

The system's wide latitude and automatic sensitivity adjustment function protects radiographs from variations in x-ray exposure conditions. This results in consistent image density for all types of radiographs. Image processing parameters are selected through anatomical region selection menus.

The FCR XL-2 reader can be network connected to multiple Flash IIP (Image & Information Processor) consoles and other FCR readers for maximum scalability and redundancy.

The FCR XC-2 supports a dedicated relationship of the reader to the Flash IIP console workstation.

A Flash IIP console workstation is provided standard with each reader unit. Fujifilm's Flash IIP is designed for simple and fast patient identification, image viewing, reprocessing and optimization. It integrates patient identification process (via manual entry or automated by connection with RIS/HIS using DICOM Worklist Management) with Fujifilm's intelligent image processing and ability to transmit images to PACS and/or imagers. The Flash IIP features a simplified user interface, which allows for complete patient processing in as few as 3 easy steps. See Flash IIP datasheet for more information.

### Options

Dedicated XL-2, XC-2 cart: Heavy duty engineered stand providing countertop and cassette holder space for the reader, with the CPU and monitor of the Flash IIP console.

**Other System Components**

- Flash IIP consoles
- Flash IIP Laptop
- FCR External DVD-RW
- Custom Stand
- Seismic Brackets
- Mobile Kit
- Axon (Sold Separately)

**Available IP Cassette Types**

- Inch settings  
14 x 17", 14 x 14", 10 x 12", 8 x 10",  
14 x 34", 14 x 40", 14 x 50", 15 x 30cm
- Metric settings  
35 x 43cm, 35 x 35cm, 24 x 30cm,  
18 x 24cm, 15 x 30cm

**Mobile kit**

Hardware brackets to stabilize & support for mobile vehicle use.  
Optional FCR XL-2 and XC-2 Laptop  
Flash IIP console also available to reduce space requirements.

**EXTERNAL DIMENSIONS AND WEIGHT\*:**

**Power Supply Conditions**

Single phase 50–60Hz  
AC 120–240V ±10%  
5A (Max.)

**Operating Conditions**

Temperature: 15 - 30°C  
Humidity: 40 - 80%RH  
(No dew condensation)

Heat Output:  
XL-995 BTU/hr  
XC-692 BTU/hr

Power consumption:  
XL-290 VA (W) or less  
XC-200 VA (W) or less

Width in. (mm)	Depth in. (mm)	Height in. (mm)	Weight lb. (kg)
23.25 (590)	15 (380)	32 (810)	218 (99)

\* Dimensions and weight are approximate and are subject to change without prior notice.

**FCR XL-2 Image Reader (Model: CR IR 359 RU2E) & FCR XC-2 Image Reader (Model: CR IR 357 RU2E)**

- (1) IP cassette type CC, LC, PC
- (2) Flash IIP Console
- (3) Image Recorder:  
DRYPIX 2000, DRYPIX 4000  
DRYPIX 5000

**Supplies**

**Fujifilm Imaging Plates (IP)**

- Standard type ST-VI  
14 x 17" (35 x 43cm), 14 x 14" (35 x 35cm), 10 x 12", 8 x 10",  
24 x 30cm, 18 x 24cm,  
15 x 30cm
- HR-V (FCR XL-2 Only)  
24 x 30cm, 18 x 24cm


**Throughput (IP/hour)**

IP/Cassette Size	XL-2	XC-2	
	Fast Scan (5pixels/mm)	Standard HQ (10pixels/mm)	Standard HQ (10pixels/mm)
14x17in	87	62	43
14x14in	94	70	50
10x12in	-	73	54
8x10in	-	87	67
24x30cm	-	75	54
18x24cm	-	92	72
15x30cm	-	73	55

**Cycle Time (Sec.)**


	XL-2			XC-2
	Standard HQ	Fast Scan	HR-V	Standard HQ
	58	41	-	82
	52	38	-	72
	49	-	-	66
	41	-	-	53
	48	-	59	66
	39	-	51	50
	49	-	-	65

**FCR XL-2/XC-2 Image Reader**

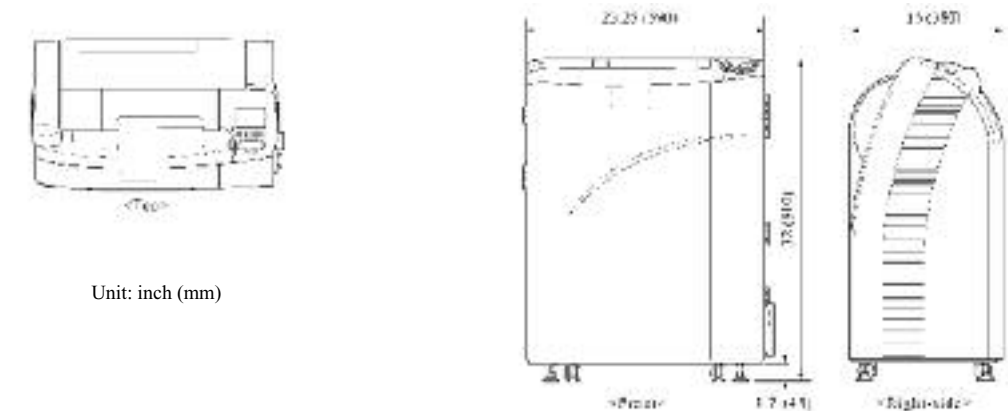

Unit: inch (mm)

**Footprint For Placement**

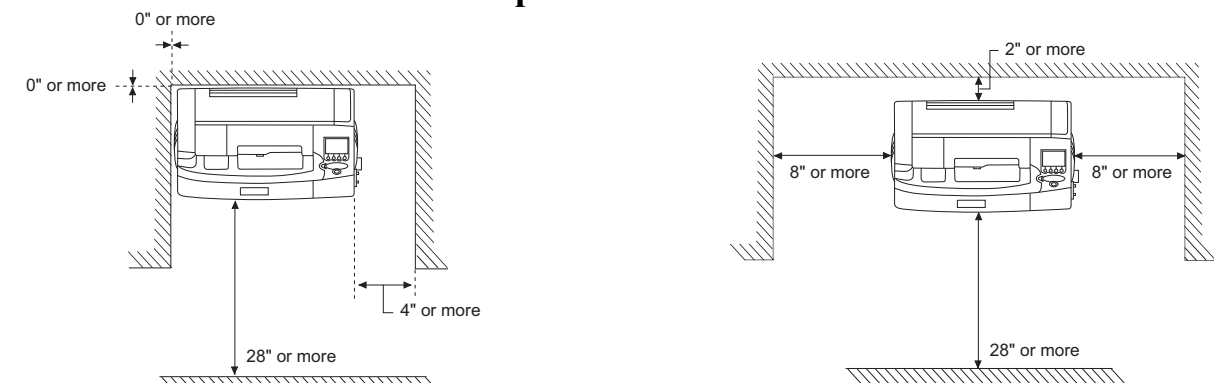

Installations Without Floor Anchoring

Seismic & Vehicle Mounted Installations  
Need 52" Minimum Height

**IMAGE READING**

Reading gray scale: 12 bits/pixel  
Output gray scale: 10 or 12 bits/pixel

**FCR XL-2 Sampling Rate**

Reading Size		14x17cm	14x14cm	10x12cm	8x10cm	24x30cm	18x24cm	HR-V 24x30cm	HR-V 18x24cm	15x30cm
Fast Scan Pixel density	Pixel Density (pixels / mm)	5	5	-	-	-	-	-	-	-
	Number of Pixels	1760x2140	1760x1760	-	-	-	-	-	-	-
High Pixel density	Pixel Density (pixels / mm)	10	10	10	10	10	10	20	20	10
	Number of Pixels	3520x4280	3520x3520	2505x3015	2000x2510	2364x2964	1770x2370	4728x5928	3540x4740	1464x2964

**FCR XC-2 Sampling Rate**

Reading Size		14x17"	14x14"	10x12"	8x10"	24x30cm	18x24cm	15x30cm
Pixel density	Pixel Density (pixels / mm)	10	10	10	10	10	10	10
	Number of Pixels	3520x4280	3520x3520	2505x3015	2000x2510	2364x2964	2364x2964	1464x2964

**Cart Dimensions**

